

Chin Implant



All implants soak in betadine and well as the pocket

Use anatomical implant Extended Anatomical Chin Implant

Place Sub periosteal


Implant should be on the lower border of the mandible

Incision should be curved

Suture implant to midline with 4 – 0 vicryl to the periosteum inferiorly

Care to avoid mental nerve – stay on the inferior border of mandible

Cheek Implant

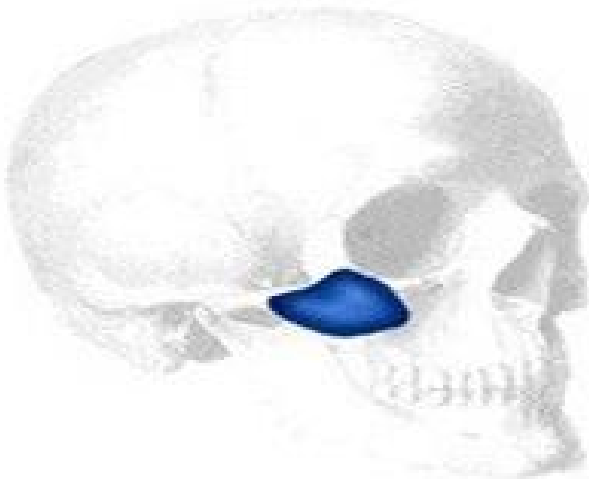

All implants soak in betadine and well as the pocket

Subperiosteal pocket

Place just above the orbital rim


Medially into the nasal malar sulcus

Fix to zygomatic buttress with a 7 mm screw

Infraorbital nerve cut a V into the medial side of implant



Temporal Implant

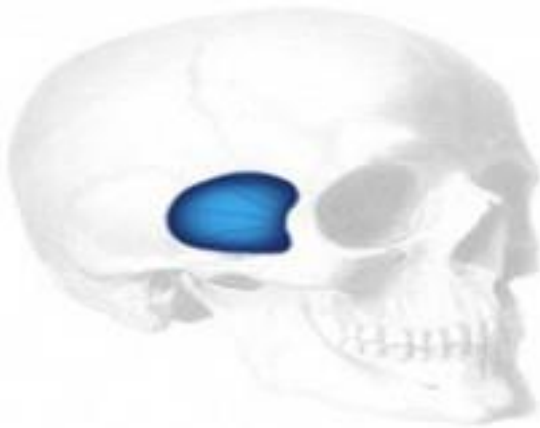

All implants soak in betadine and well as the pocket

Place under the temporalis fascia. Lays on top of temporalis muscle

Place on the temporal area , the small curve part is on the side of the orbit

Place incision a few cm behind the posterior large curve of the implant

Fix the posterior implant to the fascia to prevent displacement